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INTENDED USE
This rapid test kit is intended for the qualitative detection of SARS-CoV-2 infection from patients. It is for professional use only. It is an aid in the diagnosis of the patients with suspected SARS-CoV-2 infection in conjunction with clinical presentation and results of other laboratory tests. Results from this test kit should not be used as the sole basis for diagnosis. The test provides preliminary test results. Negative results do not preclude SARS-CoV-2 infection and should not be used as the sole basis for treatment or other management decision.

PRINCIPLE
This kit is an immunochromatography assay. According to the gold immunochromatographic test principle, double antibody sandwich method was used to detect SARS-CoV-2 antigen in the samples. When there is virus antigen presence in the sample, the antigen binds with the corresponding colloidal gold monoclonal antibody and the coated monoclonal antibody at the detection line to form a compound and then condenses into a red band, indicating a positive result. If there is no antigen in the sample, complex cannot be formed at the detection line, and no red band is shown, indicating negative result.

Whether the sample contains antigen or not, the gold monoclonal antibody will bind to the enveloped antibody at the quality control line, form a compound and condense into a red band.

- KIT COMPONENTS**
- 20 Test cassettes
 - 20 Sample tubes with prefilled sample extraction buffer
 - 20 single packaged Swabs (Nylon tip)
 - 1 Tube stand
 - 1 Instruction for use

- STORAGE AND EXPIRY**
- Store as packaged in the sealed pouch at 2-30°C, avoid hot and sunshine, dry place, valid for 24 months.
 - Do not freeze. Some protective measures should be taken in hot summer and cold winter to avoid high temperature or freeze-thaw.
 - Do not open the inner packaging until ready, it must be used in one hour if opened (Humidity ≤ 60%, Temp: 20°C-30°C).
 - Please use immediately when the humidity > 60%.

SAMPLE PRESERVATION
Saliva sample and sample of human nasopharyngeal, oropharyngeal and anterior nasal swabs should be processed as soon as possible after sample collection. If the test cannot be performed immediately, the sample should be stored in a sealed state, stored at 2-8°C for 8 hours, and stored below -20°C for 1 month. Long-term storage is not recommended.

SAMPLE COLLECTION
By anterior nasal swab:
Gently insert the entire absorbent tip of the swab (around 1.5cm) into your nostril. Firmly and slowly rotate the swab in a circular path against the inside of your nostril at least 5 times. Be sure to collect any nasal drainage that maybe present on the swab. Gently take out the swab from the nostril. Use the same swab to repeat steps in the other nostril. Then slowly take out the swab.


Saliva samples:
Do not eat, drink, smoke or chew gum for 30 minutes before collecting saliva sample. This may affect the test results. Press the tip of the tongue against the lower root of the jaw to concentrate the saliva. Place the sterile disposable swab under the tongue for at least 10 seconds. Twist it at least 5 times until the tip is completely soaked. NOTE: False negative results can occur if the saliva is not collected properly.

By nasopharyngeal swab:
The sampler gently holds the head of the person to be collected with one hand, holds the swab with the other hand, sticks the swab to the nostril to enter, and slowly penetrates backwards along the bottom of the lower nasal passage, so as not to exert too much force to avoid traumatic hemorrhage. When the tip of the swab reaches the posterior wall of the nasopharynx cavity, gently rotate it (make sure the swab contains both cells and slime), and then slowly take out the swab.


By oropharyngeal swab: The head of the person to be collected is slightly tilted and his mouth is wide open, exposing the pharyngeal tonsils on both sides. Wipe the swab across the root of the tongue.

Wipe the pharyngeal tonsils on both sides of the person to be collected back and forth with a little force for at least 3 times, and then wipe up and down the posterior pharyngeal wall for at least 3 times.

TEST STEPS
Instructions must be read entirely before taking the test. Bring test cassettes and test components at room temperature before performing the test.

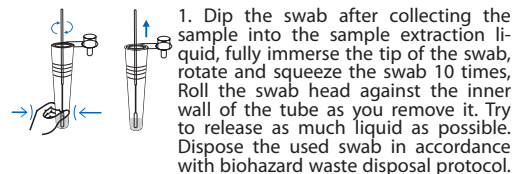2. Take out the test cassette from the sealed pouch, place it on a clean and level surface with the sample port well up. Name the test cassette.

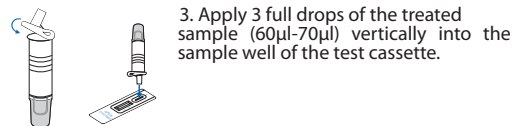4. Observe the test results immediately within 15-20 minutes, the result is invalid over 20 minutes.

Note: When using viral transport medium (VTM), it is important to ensure that the VTM containing the sample is warmed to room temperature. Cold samples will not flow correctly and can lead to erroneous or invalid results. Several minutes will be required to bring a cold sample to room temperature.

INTERPRETATION OF RESULTS

POSITIVE:
Two (2) distinct colored lines appear. One line should be in the control region (C) and the other line should be in the test region (T).

NEGATIVE:
One (1) colored line appears in the control region (C). No apparent colored line appears in the test region (T). The negative result does not indicate the absence of analytes in the sample.

INVALID:
No colored lines appear, or control line fails to appear, indicating that the operator error or reagent failure. Verify the test procedure and repeat the test.

- LIMITATIONS**
1. This reagent is a qualitative detection reagent, which cannot determine the exact content of antigen.
 2. The test results of this reagent are only for the reference of clinicians and should not be taken as the sole basis for clinical diagnosis and treatment. Clinical management of patients should be considered in the light of their symptoms/signs, medical history, other laboratory tests and treatment responses.
 3. Restricted by antigen detection reagent method, the lowest detection limit (sensitivity analysis) is generally lower than that of nucleic acid detection, so the researchers deal with negative result to give more attention, should be combined with other test results comprehensive judgment, advice to doubt the negative result of nucleic acid detection or virus isolation culture identification method for review.
 4. False negative results may be caused by unreasonable sample collection, transport and treatment, and low viral load in samples.

PERFORMANCE CHARACTERISTICS

1. Sensitivity and Specificity
The verification of samples of all collection methods is carried out with RT-PCR and the results meet the following requirements:
The National Reference Panel for SARS-Cov-2 Antigen Detection Kit of the National Institutes for Food and Drug Control was used for verification, and the results met the following requirements:

Coincidence rate of negative reference products: National negative reference materials (N1-N20) were tested, and the results were all negative, with a coincidence rate of 100%.
Coincidence rate of positive reference products: tested with national positive reference products (P1-P8), the results are all positive, and the coincidence rate is 100%.

Repeatability: The 10 test results of R1 and R2 should be positive, and the color rendering is uniform.

Clinical Sensitivity and Specificity of naso-oropharyngeal samples:

SARS-CoV-2 Antigen Rapid Test		RT-PCR		Total
		Positive	Negative	
SARS-CoV-2 Antigen Rapid Test	Positive	243	9	252
	Negative	7	541	548
Total		250	550	800

Diagnostic sensitivity: 97.2% (94.3% - 98.9%)*
Diagnostic specificity: 98.4% (96.9% - 99.2%)*
Total coincidence rate: 98.0% (96.8% - 98.8%)*
*95% Confidence Interval

Clinical Sensitivity and Specificity of anterior nasal swab samples:

SARS-CoV-2 Antigen Rapid Test		RT-PCR		Total
		Positive	Negative	
SARS-CoV-2 Antigen Rapid Test	Positive	146	1	147
	Negative	4	149	153
Total		150	150	300

Sensitivity: 97.33% (93.31%-99.27%)*
Specificity: 99.33% (96.34%-99.98%)*
Total consistent: 98.33% (96.15%-99.46%)*
*95% Confidence Interval

Clinical Sensitivity and Specificity of saliva samples:

SARS-CoV-2 Antigen Rapid Test		RT-PCR		Total
		Positive	Negative	
SARS-CoV-2 Antigen Rapid Test	Positive	95	1	96
	Negative	5	99	104
Total		100	100	200

Sensitivity: 95.0% (88.72%-98.36%)*
Specificity: 99.0% (94.55%-99.97%)*
Total consistent: 97.0% (93.58%-98.89%)*
*95% Confidence Interval

2. Limit of detection:

LOD concentration	30 TCID ₅₀ /mL
-------------------	---------------------------

3. Cross-reactivity
With human coronavirus 229E, human coronavirus OC43, human coronavirus NL63, adenovirus, human metapneumovirus, parainfluenza virus 1, parainfluenza virus 2, parainfluenza virus 3, parainfluenza virus 4, influenza A, type B influenza, enterovirus, respiratory syncytial virus, rhinovirus, SARS coronavirus, MERS coronavirus, Haemophilus influenzae, Streptococcus pneumoniae, Streptococcus pyogenes, Candida albicans, Bordetella pertussis, Mycoplasma pneumoniae, pneumonia Chlamydia, Legionella pneumophila, etc. have no cross reaction (see table in the German IFU).

4. Interfering
Common interfering substances in the sample, such as blood, mucin, and pus, have no effect on the test results.

- PRECAUTIONS**
1. For IN VITRO diagnostic use only.
 2. Reagents should be used as soon as possible after opened. This reagent cannot be reused for disposable.
 3. The test device should remain in the sealed pouches until use.
 4. All specimens and reagents should be considered potentially hazardous and handled in the same manner as an infectious agent after use.
 5. Don't use after the expiration date.

- LITERATURHINWEISE | REFERENCES**
1. Forni, D., Cagliani, R., Clerici, M. & Sironi, M. Molecular evolution of human coronavirus genomes. Trends Microbiol. 25, 35-48 (2017).
 2. Ithete, N. L. et al. Close relative of human Middle East respiratory syndrome coronavirus in bat, South Africa. Emerg. Infect. Dis. 19, 1697-1699 (2013).

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